



EquiSys – Safety and Style

EquiSys was designed for the safe mounting of all technological devices within the medical environment but with a style that enhances, rather than detracts. Rated for a safety factor of 6, you can be confident that your installation is safe and secure.

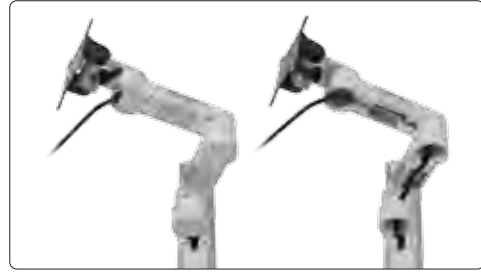

No other arm system currently on the market features such complete cable integration as the EquiSys line, further enhancing safety and aesthetics (shown with covers in place, and removed).

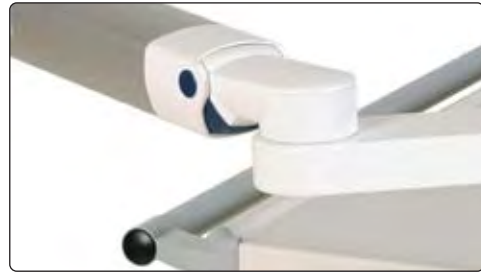

When designing the EquiSys product line, great emphasis was placed on the fit and finish for cleaning as well as aesthetic reasons. Smooth, powder coated anti-microbial finishes ensure that proper disinfection procedures can be carried out both easily and effectively.

EquiSys Sample Configuration #1

This configuration features a twin flat panel assembly allowing each panel to tilt and pivot independent of one another. The double swing arm keyboard tray allows the keyboard to be drawn straight out and then tucked neatly under the flat panels when not in use. The power supply holster completes the system. Potential application: PACS.

- #65-D910G-4 Wall mount, dual flat panel VESA 75/100
- #65-24240E-4 Wall mount, keyboard tray on double swing arm
- #65-14-4 Wall mount, universal power supply holster
- #65-0096 Wall channel, 960mm length with covers



EquiSys Sample Configuration #2

This system features a wall mounted variable height adjustable arm with extension arm for optimal height adjustment whether standing or sitting. The small wire basket can be used as storage. Potential application: Patient bedside data entry.

- #65-53024F-4 Wall mount, variable height arm and extension arm VESA 75/100
- #65-52-4 Wire basket, small 200mm x 150mm x 100mm (W x D x H)
- #65-0096 Wall channel, 960mm length with covers



EquiSys Sample Configuration #3

The upper arm features a variable height arm with extension arm providing exceptional reach from the mounting point as well as height adjustment for optimal viewing. The lower single swing arm allows a second flat panel to be viewed and to pass underneath the upper arm. Potential application: Laparoscopy and video carts.

- #65-53024G-4 Wall mount, variable height arm and extension arm VESA 75/100
- #65-230G-4 Wall mount, single swing arm
- #65-0020 Wall channel, 200mm length with covers

EquiSys Sample Configuration #4

This configuration features a flush, flat panel mount which provides tilt and pivot of the flat panel. A power supply holster is shown mounted. Also featured is our compact Medi-Key keyboard. Wireless Medi-Key can be disengaged from its mount for use as well as cleaning. Potential application: Data entry.

- #65-910G-4 Wall mount, flush tilt/pivot flat panel mount VESA 75/100
- #65-950E-99 Fold up keyboard bracket for Medi-Key keyboard
- #65-14-4 Wall mount, universal power supply holster
- #65-0096 Wall channel, 960mm length with covers
- #SMK-A41-01-US Wireless illuminated keyboard with integrated touch pad, disinfectable




Keyboard shown in folded "up" position.

EquiSys Sample Configuration #5

Ideal for small device mounting, this configuration shows a single swing arm with universal device mounting tray. Quick release as well as tilt and pivot of tray are standard. Large wire basket with integrated cable hooks provides useful storage capabilities. Potential application: Patient monitoring.

- #65-224A-4 Wall mount, single swing arm with universal mounting plate
- #65-51-4 Wire basket, large 300mm x 200mm x 180mm (W x D x H)
- #65-0072 Wall channel, 720mm length with covers

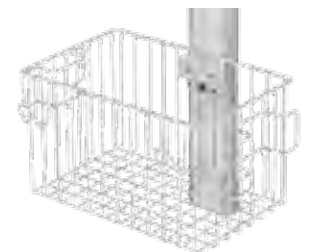




Surgimed Solutions, LLC

EquiSys
Equipment Mounting System

EquiSys Equipment Mounting System

EquiSys is modular in design allowing for literally hundreds of configurations. A few examples have been shown in order to spur thought and to provide a basis from which to begin. For further information, please visit our website or contact us at the number below in order to discuss this or any of our other fine products.

The EquiSys product is designed and manufactured by CIM Med GmbH of Germany. Surgimed Solutions, LLC is the US Distributor for CIM Med GmbH.

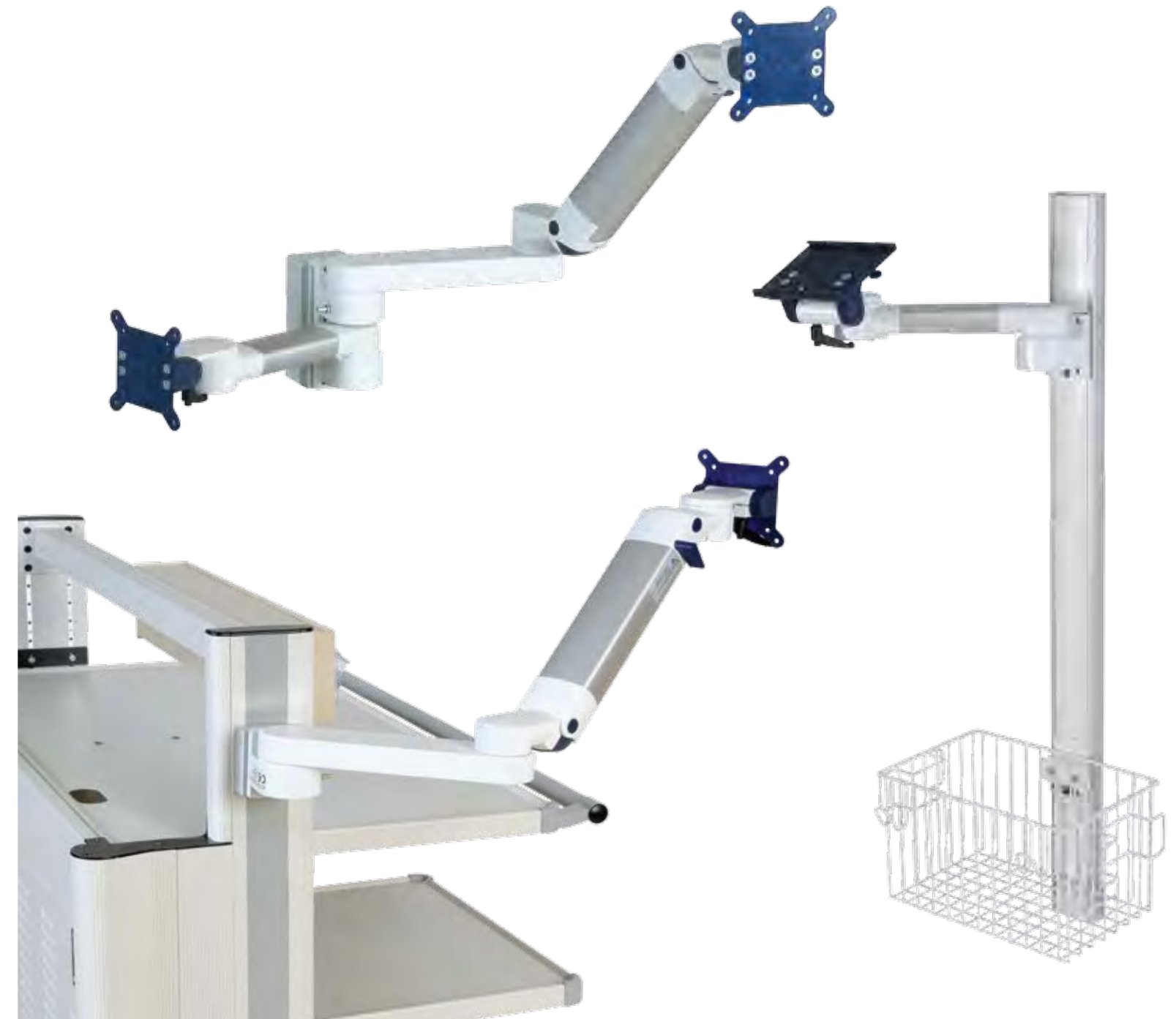

Legacy
Carts & Roll Stands



P E N T A
Medical Grade PCs



Medi-Key™
Disinfectable Keyboards



Surgimed Solutions, LLC